

Soredex Cranex D

Integration Instructions for ViewPoint



The Soredex Cranex D panoramic dental X-ray machine, when set up correctly, will function the same as a regular scanner within VP Imaging. No additional hardware or software is needed.


Setup

Install the Soredex Cranex D, appropriate software, and corresponding TWAIN driver. This should be accomplished by someone from Soredex.

For information on setting up VP Imaging, consult the VP Imaging documentation and contact Ortho2 Software Support with any questions.

Using the Integration

- 1 Open a patient record in ViewPoint and click the **VP Imaging** tab.
- 2 Select the appropriate layout.
- 3 Click **Scan** to open the *VP Imaging - Scan* window.
- 4 Select the appropriate timepoint.
- 5 Select the appropriate image type.
- 6 Click **Select Scanner** to open the *Select Scanner* window.
- 7 On the **Twain** tab, select the Soredex device. Click **OK**.
- 8 Click **Scan** on the *VP Imaging - Scan* window. This will start your scanning software. For more information on proceeding, refer to the Soredex documentation.
- 9 A **Count=1** notation will appear beneath the **Close** button. Scan additional images or click **Close**.

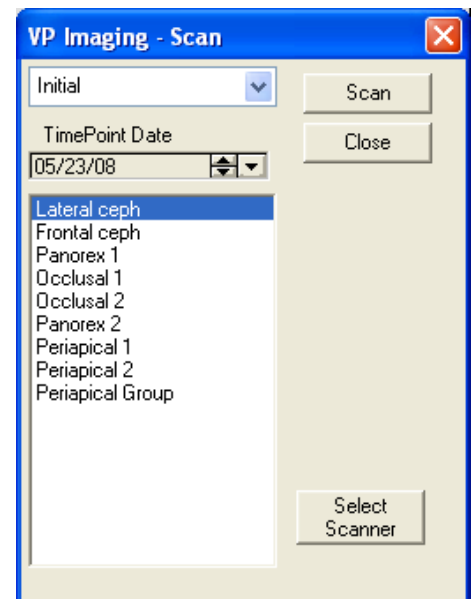

Refer to the VP Imaging documentation for instructions on how to use the Image Editor.

Contact

Soredex Cranex D
Soredex
300 West Edgerton Avenue
Milwaukee, WI 53207
www.soredexusa.com
(800) 235-8854

ViewPoint
Ortho2
1107 Buckeye Avenue
Ames, IA 50010
www.ortho2.com
(800) 346-4504 or (515) 233-1026