

# Visix

Integration Instructions for ViewPoint




The ViewPoint integration with Visix software from Air Techniques provides direct access to the Visix program from ViewPoint. If the patient doesn't exist in Visix, one will be created with first name, last name, and patient ID from the patient's ViewPoint record. These features require the ViewPoint Third Party Integration.


## Setup

Before you begin, be sure both ViewPoint and Visix are installed properly.

- 1 From the ViewPoint Main Menu, click **[Tools & Utilities]** and then **[Third Party Integration]**.  
*Note: If you receive a CD-Key mismatch message, contact your Ortho2 Systems Consultant.*
- 2 Click **[New]**.
- 3 Select **Command Line** for the integration type.
- 4 Enter "Visix" as the **Integration Name**.
- 5 For the **Filename**, use [...] to navigate to the **Tiger1.exe** and select it.
- 6 Enter "\*.VISIX" (all caps) as the **Command Line Suffix**.
- 7 Click **[Add]** and **[OK]**.

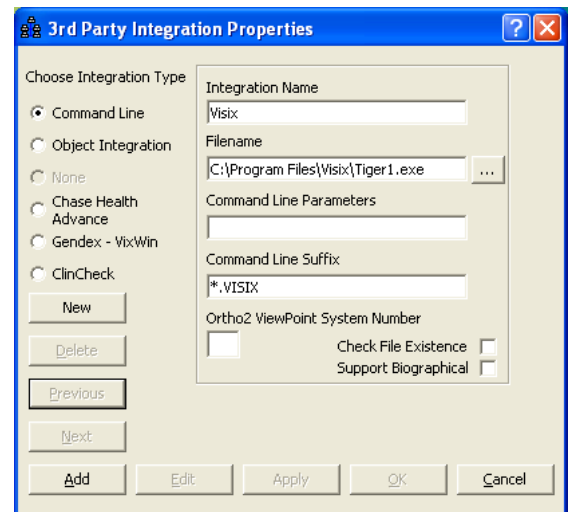

## Using the Integration

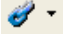To open Visix from a ViewPoint patient record, simply click the **Third Party Integration**  button in the ViewPoint toolbar and select **Visix**. Visix will start and automatically load the current patient. If no record for this patient currently exists, one will be created with the first name, last name, and patient ID.

With this integration, the 3rd Party **Update/New** button  in the patient folder has no function.

## Contact

### Visix

Air Techniques, Inc.  
1295 Walt Whitman Road  
Melville, NY 11747  
www.airtechniques.com/home.aspx  
800 AIR-TECH

### ViewPoint

Ortho2  
1107 Buckeye Avenue  
Ames, IA 50010  
www.ortho2.com  
(800) 346-4504 or (515) 233-1026